

MEGA NEWSTM Hcg Pregnancy Test








50 Tests
Pack



WELCOME TO CLARITY


Discover the accuracy and simplicity of the **Mega News HCG Pregnancy Test Kit**, your reliable partner in early pregnancy detection. With our advanced technology, you can achieve peace of mind with just a few simple steps.

Why Choose MEGA NEWSTM ?


 <h3>Early Detection</h3> <p>Detect pregnancy hormones as early as 5 days before your missed period.</p>	 <h3>High Accuracy</h3> <p>Over 99% accuracy in detecting pregnancy when used from the first day of your missed period.</p>	 <h3>Ease of Use</h3> <p>Simple and user-friendly for quick, at-home testing.</p>	 <h3>Fast Results</h3> <p>Clear and definitive results in just 3 minutes.</p>	 <h3>Discreet & Private</h3> <p>Perform the test in the comfort of your home, with no need for doctor visits or appointments.</p>
---	--	--	---	--

How it work

Test Procedure

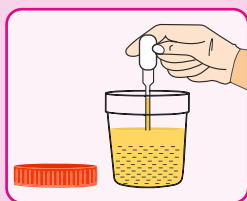


1. Collect Sample:
Using the provided dropper, collect a sample of your urine.

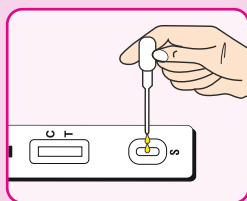

2. Apply Sample:
Add 2 drops of urine to the designated area on the test strip.




3. Wait:
Allow the test to process for 5-10 minutes.

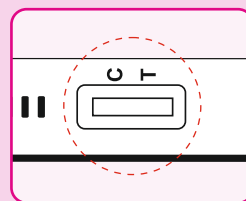


4. Read Results:
Check the result window for clear indicators of pregnancy.


Interpretation of results

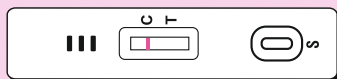

NEGATIVE:
One pink line appears.
The test indicates no pregnancy.

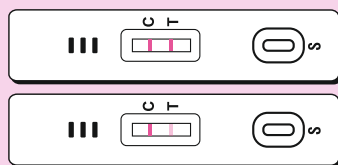

POSITIVE:
Two pink lines appear.
The test confirms pregnancy.

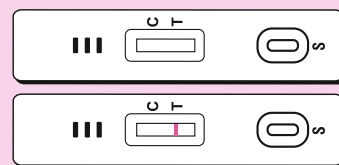

INVALID:
No lines or an unclear result may indicate an error.
Please read instructions carefully and try again if necessary.

Frequently Asked Questions

• How accurate is the test?

Mega News HCG Pregnancy Test Kit offers over 99% accuracy when used as directed.

• When should I take the test?

For the most accurate results, use the test after your missed period. Early testing can be done 5 days before a missed period, but accuracy improves with time.

• Can I use the test with any urine sample?

For best results, use the first-morning urine, which has the highest concentration of the pregnancy hormone.

• What should I do if the test result is positive?

Follow up with a healthcare provider for confirmation and to discuss next steps.

Product Details

• Contents: 1 individually wrapped test card, 1 dropper • Shelf Life: 24 months from the manufacturing date • Storage: Store in a cool, dry place

Mega News HCG Pregnancy Test Kit – Your trusted solution for early and accurate pregnancy detection. For more information, visit our website or contact our support team.

Mega News – Empowering you with clarity and confidence.

Manufactured by:

OPTSARASMED (OSD) HEALTHCARE PVT. LTD.

(ISO 13485:2016 Certified Company)

H1-420, Khushkhera, RIICO Industrial Area, Alwar, Bhiwadi-301707, Rajasthan, India

Customer Care No.: 01493-403808 | E-mail: marketing@osdhealthcare.in

Website: www.osdhealthcare.in